„Airway Alex“ Wendl/Guedel-Trainer

Merkblatt / Factsheet

Art. Nr. / Item No. 7930

Der Airway Trainer „Airway Alex“

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
The “Airway Alex” airway trainer

has been specifically designed for hands-on airway management training. This compact training model allows realistic practice in the insertion of nasopharyngeal and oropharyngeal airways and is ideal for the education and training of EMS professionals, tactical medicine, nurses, medical staff and first aid providers.

The flexible face features anatomically designed nasal openings for the insertion of a nasopharyngeal airway (NPA). The trainer provides two different anatomical conditions: The left nostril has a normal opening, while the right nostril is narrower and simulates a deviated nasal septum. This allows users to practice different levels of difficulty when inserting a nasopharyngeal airway.

The open mouth allows both visual confirmation of the NPA position and the insertion of an oropharyngeal airway (OPA). The realistically reproduced airway extends to the epiglottis and provides realistic tactile feedback during insertion. A soft, movable tongue further enhances the training experience and simulates patient anatomy.

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
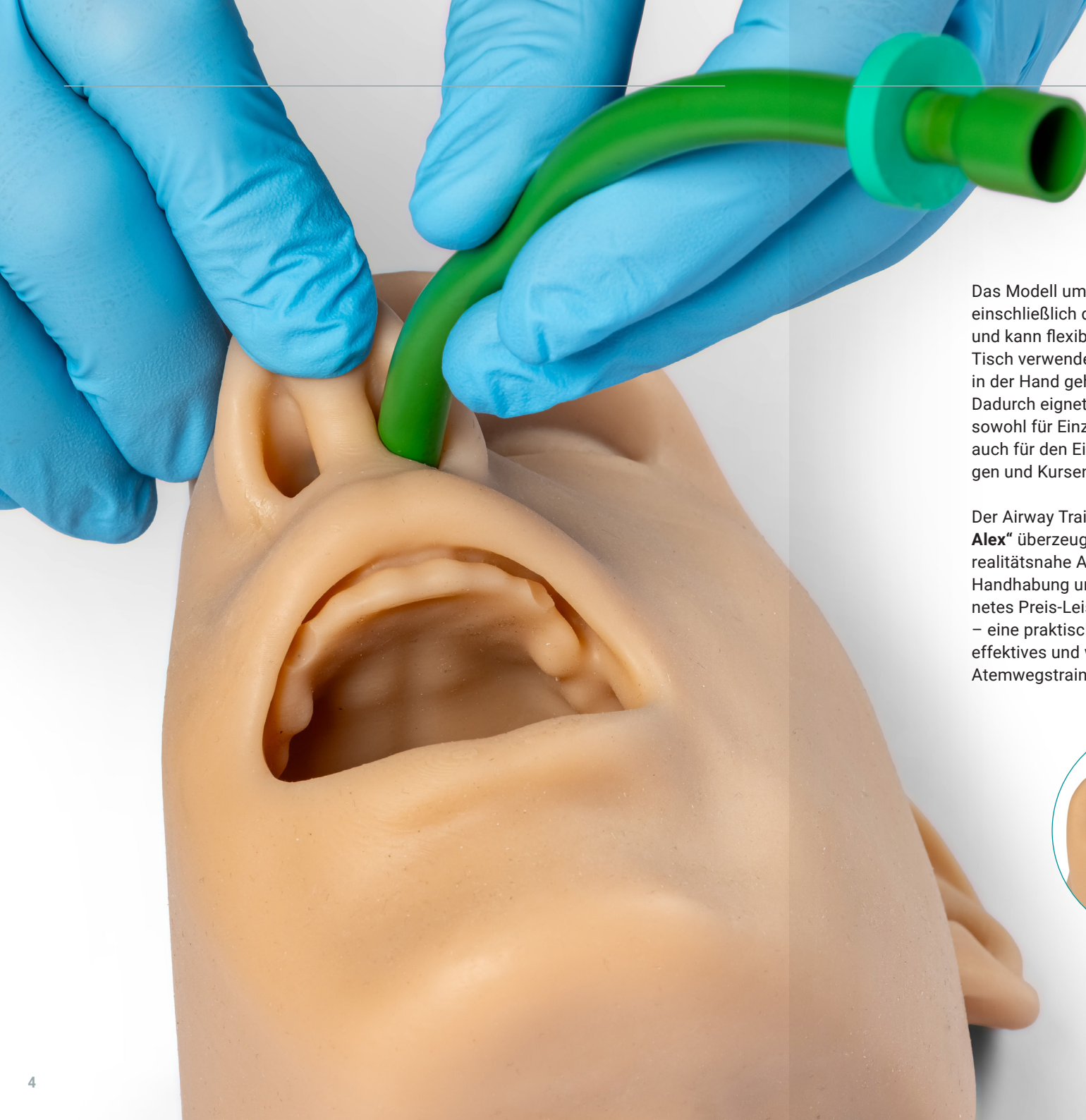
The model consists of the face including the left ear and can be used conveniently on a table or held in the hand. This makes the trainer suitable for both individual practice and use in professional training courses.

*The “**Airway Alex**” Airway Trainer combines realistic anatomy, easy handling and excellent value for money – a practical solution for effective and repeatable airway management training.*

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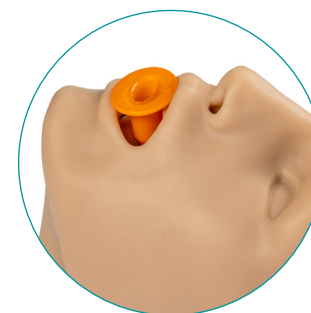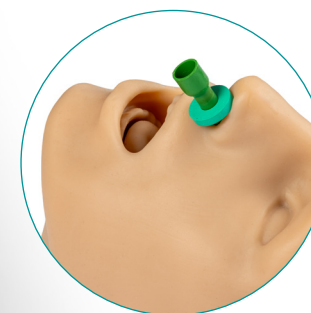Product Features

- Realistic airway trainer for airway management
- Anatomically accurate airway extending to the epiglottis
- Different nasal anatomies for varying levels of difficulty
- Narrower right nostril simulating a deviated nasal septum
- Left nostril with normal opening
- Flexible facial structure with a soft, movable tongue
- Nasal openings for inserting a Wendl tube
- Open mouth for checking the position of the Wendl tube and for inserting a Guedel airway
- Face includes the left ear
- Suitable for tabletop or handheld exercises
- Ideal for emergency medical services, emergency medicine, nursing, training, and simulation
- Cost-effective training model with high practical relevance

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Ideal for training in:

- Insertion of a nasopharyngeal tube
- Management of different anatomical conditions during nasopharyngeal airway insertion
- Positioning of the nasopharyngeal tube
- Insertion of a Guedel airway (oropharyngeal tube)
- Anatomical orientation in the upper airway
- Basics of airway management

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Product Features:

After use, clean silicone surfaces with a soft, damp cloth and a mild detergent. Then let them dry completely. Do not use harsh cleaning agents, disinfectants, solvents, abrasive cleaners, or alcohol-based products, as these can damage the material.

Silicone is durable, but not indestructible. Avoid contact with sharp-edged or pointed objects such as watches, jewelry, or similar items, as these can cause damage.

Store the product in a dry, clean place, protected from direct sunlight.

Use only fake blood that does not permanently stain the skin.

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Erler-Zimmer Medical GmbH